

Wahl® Heated Lumbar Massage™ Cushion Quickstart Guide



1
Insert adaptor jack into the plug on the right side of the unit.



2
Plug wall adaptor into AC 120V outlet or car adaptor into cigarette lighter power supply.



3
To turn the unit on, press the power button on the controller once.



4
To turn the heat on, press the heat button once. Press again to turn heat off.



5
To change the programs, press the power button. Select from multiple settings by pressing the programs button successively until desired program is found.



6
To turn the unit off, hold the power button for 3 seconds. This unit automatically turns off after 20 minutes.


Wahl® Heated Lumbar Massage™ Cushion Key Features

Prepare for use at home, in the office, or in the car:

Place the unit in a chair or car seat, and pull the attached elastic strap over top of the chair or seat to hold the cushion in place. Make sure the unit is flush with the chair or seat to allow for maximum benefit of the lumbar support.

Programs

LOW Vibration — for soothing massage
HIGH Vibration — for vigorous, deeper massage
OFF — if you prefer heat only

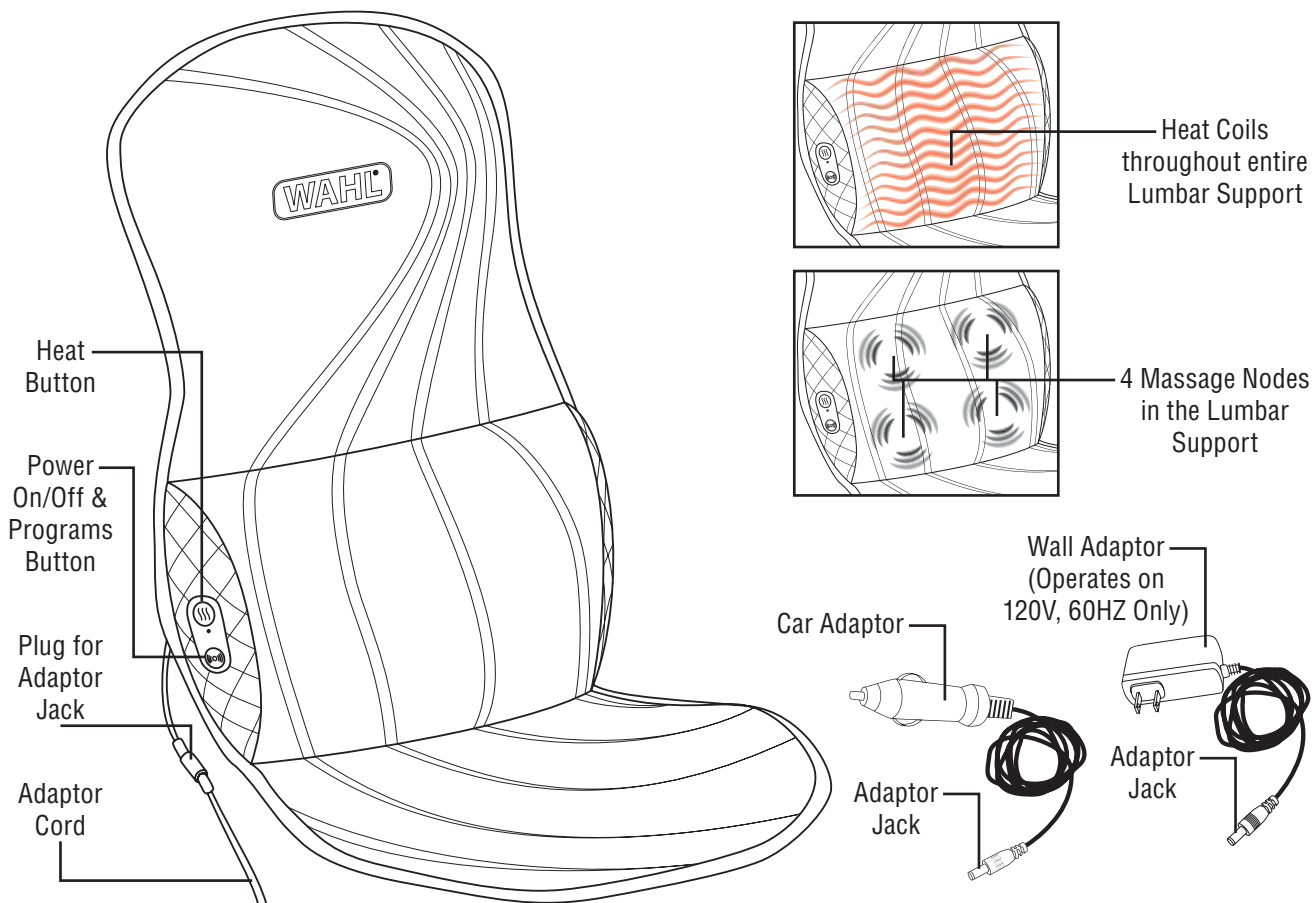

If you are not completely satisfied with your Wahl Massager, please contact us before returning this product at WahlUSA.com or 1-800-767-9245.