

# Wahl® Vibrating Heat Neck & Back Wrap Quickstart Guide



Insert adapter jack into the plug on the left side of the product.



Plug into a AC 120-volt outlet.



To turn the unit on, press the power button on the controller once. Press again to switch between high speed and low speed.



To turn the heat on, press the heat button once. Press again to turn it off. Heat may be used alone or with massage.




To change the programs, press the program button. Select from multiple settings by pressing the programs button successively until desired program is found.



To turn the product off, hold the power button for 3 seconds. This unit automatically turns off after 20 minutes.

# Wahl® Vibrating Heat Neck & Back Wrap Key Features

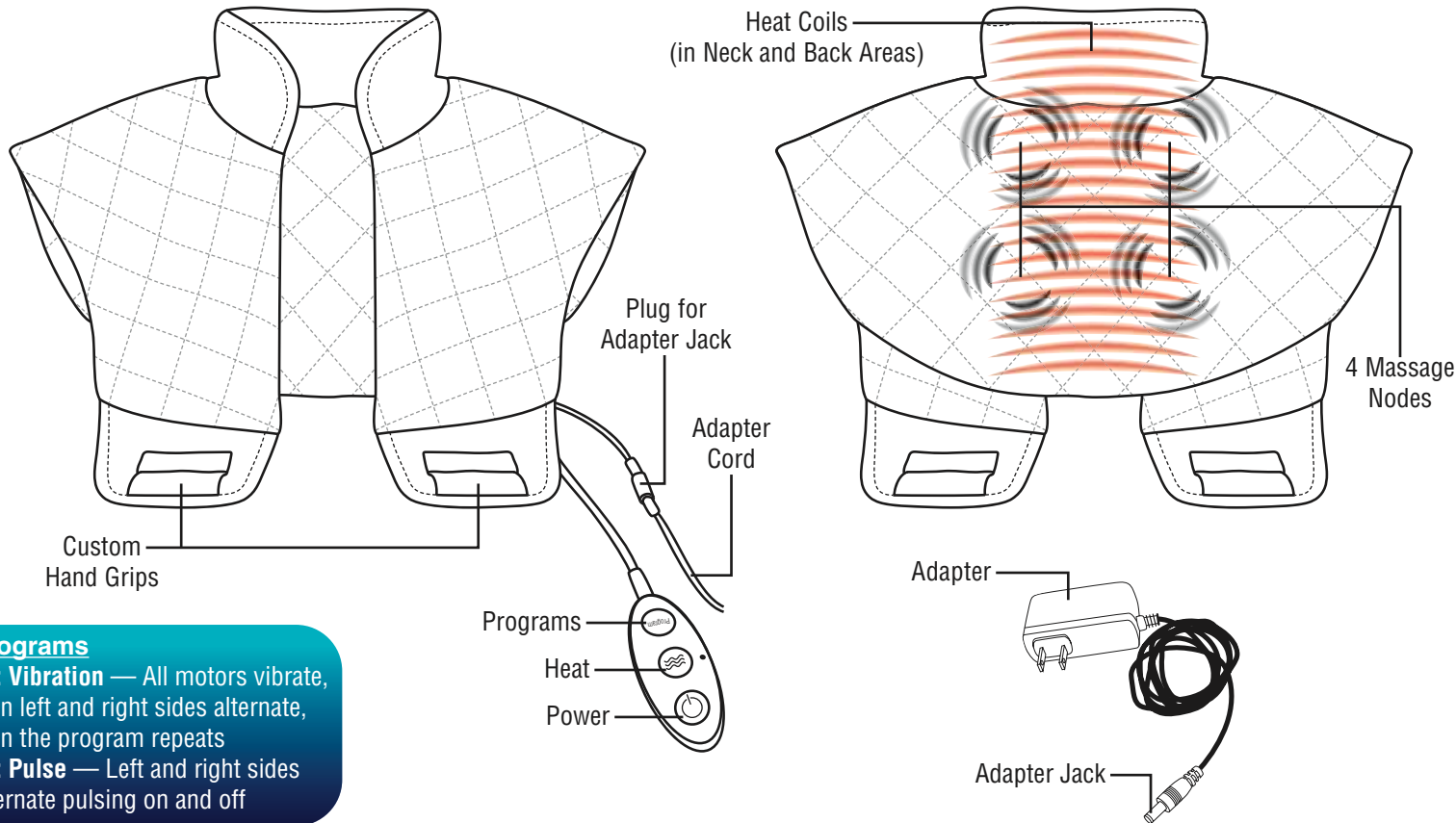

## Programs

**#1: Vibration** — All motors vibrate, then left and right sides alternate, then the program repeats

**#2: Pulse** — Left and right sides alternate pulsing on and off